



Urinalysis

from Siemens Healthineers



reddot winner 2025




TIMIK present the total solution for urinalysis from SIEMENS Healthineers in the Nordics. The Atellica LumIQ Analyzer presents next generation of Urinalysis, combining the proven precision of the CLINITEC family with advanced digital imaging technology.

Trusted chemistry. Advanced technology. Modern design. Illuminating the intelligence of urinalysis.

Atellica LumiQ Analyzer is a portable, easy-to-use, semi-automated analyzer designed to read Siemens Healthineers Multistix and hCG cassettes, producing results intended for in vitro diagnostic use only in point-of-care and clinical laboratory settings. Its compact, human-centric design makes it easy to move and use, while the modern, intuitive interface and digitised workflows reduce the likelihood of user error.



Test table turn for hCG or stix analysis

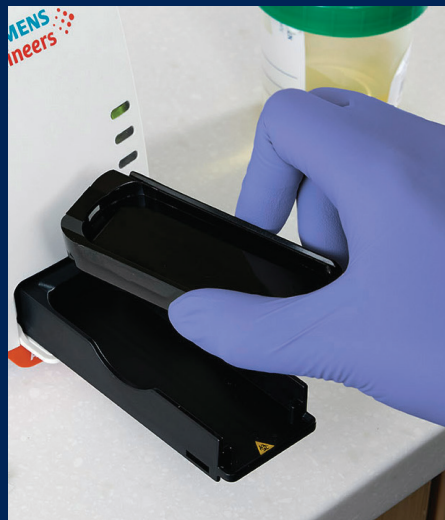

Atellica LumiQ is a portable and easy-to-use urine analyser that performs both:

- Urine strip testing (Multistix)
- hCG pregnancy testing (cassette)

Scalability

Atellica LumIQ integrates seamlessly into existing infrastructures.

The analyser's unique three-to-one module-to-display ratio allows for simultaneous testing, making it easier to scale capacity and meet increasing demands while reducing total costs of ownership.



reddot winner 2025



Enhanced UI

Modern, intuitive design improves the user experience.




Digital safeguards

Integrated camera and advanced algorithm, reduce human subjectivity and minimise user errors.



Reliable results

Routine urinalysis, microalbumin, and hCG pregnancy testing.

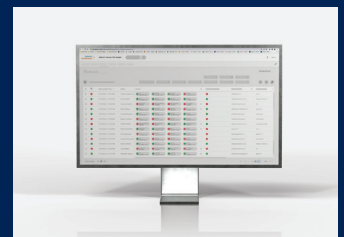

Built-in connectivity

Secure link to LIS/HIS/EHR via ethernet, Wi-Fi, or middleware.



Effortless QC

Convenient, liquid-ready Chek-Stix UA Dipeptide Quality Control.



Compact analyzer

Smaller footprint when space is at a minimum.



Your workflow, your way

Add more testing modules when the need arises.



Portability


Plug-and-play functionality to bring results directly to the patient.

Well known Multistix

– in a brand-new look!

Auto-Checks technology

Built in Auto Checks technology reduces subjectivity and human error by providing consistent, reliable insights in seconds.

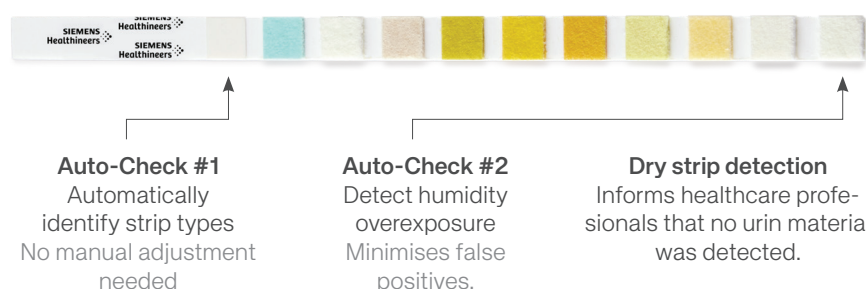

Specifications Atellica LumIQ

Time to prepare test:
12 seconds

Time to result:
45 seconds

Test menu:
Automatic read of Siemens Healthineers Multistix 5, 7, 8 SG, 10 SG, Uristix, Albustix, microalbumin and CLINITEST hCG.

Calibration:

- Calibrated during manufacturing
- The test table's calibration bar is read before each test

Connectivity:

- Build in Wi-Fi
- Lan cable via adaptor accessory

Data storage:

- QC results: 5,000 results
- Patient results: 10,000 results
- Operator list: 10,000 operators

Environmental requirement Atellica LumIQ:

- 15°C to 32°C
- 10% to 90% RH (noncondensing)
- 4500 m (14,764 ft)
- Keep away from direct sunlight



Dipeptide Liquid QC

3. Party quality control material, distributed by Siemens Healthineers

